

Whole slide and microscopy image analysis with QuPath and OMERO

Léo Leplat¹, Alan O'Callaghan¹, and Peter Bankhead¹

¹ Institute of Genetics and Cancer, University of Edinburgh, Edinburgh, EH4 2XU, UK
Corresponding author

DOI: [10.xxxxxx/draft](https://doi.org/10.xxxxxx/draft)

Software

- [Review](#)
- [Repository](#)
- [Archive](#)

Editor: [Erick Martins Ratamero](#)



Submitted: 16 March 2026

Published: unpublished

License

Authors of papers retain copyright and release the work under a Creative Commons Attribution 4.0 International License ([CC BY 4.0](#))

Summary

QuPath is open-source software for bioimage analysis. As a desktop application that is flexible and easy to install, QuPath is used by labs worldwide to visualise and analyse large and complex images. However, relying only on images stored only on a local file system limits QuPath's use for larger studies.

This paper describes a new extension that enables QuPath to access pixels and metadata from an OMERO server. This enhances the software by allowing it to work efficiently with images stored remotely, while also serving as a template for developers who want to connect QuPath to other image management systems.

Statement of need

QuPath is a popular, open-source application for image analysis, written in Java (Bankhead et al., 2017). With over 1,000,000 software downloads and more than 6,000 paper citations to date, QuPath's support for large and complex images have helped establish the software as a key biomedical research tool.

QuPath is routinely used to analyse whole slide images, which are common in research and the field of digital pathology. A 'small' whole slide image might be 120,000 × 60,000 pixels in size, which equates to around 20 GB uncompressed data (assuming 8-bit, 3-channel RGB pixels). Some images used with QuPath can be much larger, such as fluorescence multiplexed whole slide images that often contain dozens of 16-bit or 32-bit channels. While lossy compression can somewhat reduce file sizes, data management remains a major issue for large studies.

OMERO is a popular, open-source image management solution that enables images to be stored on a central server and viewed through a web browser (Allan et al., 2012). It is installed in institutions worldwide and is a key technology powering the Image Data Resource (IDR) (Williams et al., 2017) — a major repository currently hosting over 400 TB of published imaging data. OMERO also supports whole slide and multiplexed images, making it a natural fit for many QuPath use cases.

The QuPath OMERO extension bridges the gap between both tools, making it possible to apply QuPath analysis to images hosted in OMERO. By efficiently accessing only the required pixels and metadata, the extension avoids the need to download and duplicate entire datasets.

State of the field

The current work builds upon lessons learned in the development of two previous QuPath extensions that connected to OMERO.

Initially, within the QuPath development team we created the [QuPath Web OMERO extension](#). This lightweight extension could be added to QuPath and used to access images through the OMERO web API. As a single jar file without external dependencies, this was easy to install and use. However, it had the major limitation of being able to request only JPEG-compressed, 8-bit RGB tiles — not raw pixel data. This made it unsuitable for quantitative analysis and incapable of supporting most fluorescence images, which can have different bit-depths and channels.

The [BIOP OMERO extension](#) started as a fork of our original extension. It was developed by the [BioImaging and Optics Platform \(BIOP\)](#) team to use the OMERO Ice API, rather than the web API. This had the major advantage of raw pixel access, but also had disadvantages: it required many additional dependencies, raw pixel access could be slow, and it supported only OMERO servers with authentication.

We developed the present QuPath OMERO extension from scratch to provide a more flexible and maintainable solution. It aims to overcome previous limitations, while adding support for recent OMERO features and implementing best practice. Key features include the ability to:

1. Browse and import images from both public and private OMERO servers.
2. Access pixels flexibly using different APIs.
3. Exchange extra information (e.g., annotated regions of interest) between QuPath and OMERO.
4. Run custom scripts within QuPath to interact with the OMERO server.

The new extension is also the first to include unit tests. It can be easily installed through QuPath's new extension manager, which we developed in parallel to support downloading both the extension and its optional dependencies.

Software design

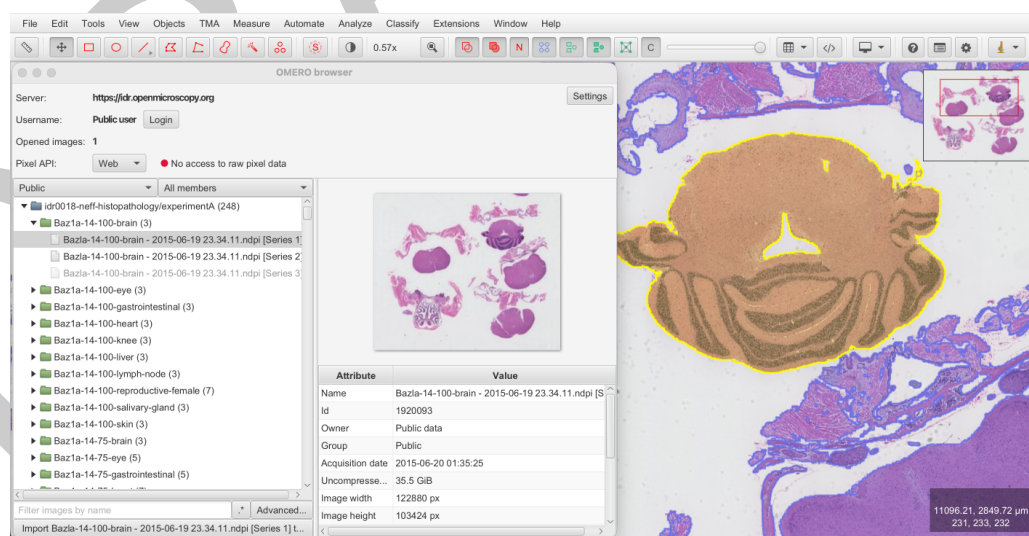


Figure 1: Screenshot showing the QuPath OMERO extension browsing the [IDR](#).

While the need for flexible pixel access from both public and private OMERO servers was the initial motivation for this work, we developed a new extension (rather than another fork) to avoid the constraints of earlier design choices and backwards compatibility. This allowed us to better separate core logic from the user interface, improve responsiveness, and add extensive tests. This separation is reflected in the code being divided into two main packages: `gui` and `core`.

68 The `gui` package contains the user interface controls and interaction with the main QuPath
69 application. It is designed to follow the general style of QuPath and to always stay responsive
70 by delegating long-running tasks to background threads. When browsing an OMERO server,
71 the interface resembles the OMERO web client to streamline user adoption.

72 The core package is largely independent of the user interface, although we permitted the use
73 of JavaFX observable properties and some QuPath classes (e.g., to manage preferences). This
74 was a pragmatic decision that allowed us to write simpler code with less duplication. The
75 key distinction between packages is therefore that core contains classes that can be used
76 headlessly and via scripts, while `gui` requires the interface to be instantiated.

77 The core package is also the focus of our unit testing. Because many functions can only be
78 tested when an active connection to an OMERO server is established, a Docker container that
79 hosts an OMERO server is automatically created when unit tests are run. The extension also
80 provides a bash script to create this Docker container outside of the unit tests, which can be
81 useful for manual testing.

82 A central part of the core package is support for different pixel APIs. Three are currently
83 implemented:

- 84 ■ The **web** pixel API: this is enabled by default and is available on every OMERO server.
85 It is fast, but suffers from the same limitations as the QuPath Web OMERO extension:
86 only RGB images can be read, and images are JPEG-compressed.
- 87 ■ The **Ice** pixel API: similar to the BIOP OMERO extension, this can read any image and
88 access raw pixel values. It is available when its additional dependencies are installed and
89 when the server requires authentication.
- 90 ■ The **pixel data microservice** API: this can read any image and access the raw pixel values.
91 It works for both public and private servers. It requires no additional installation on the
92 client side, but requires that the [OMERO Pixel Data Microservice](#) is installed on the
93 server.

94 The variety of server configurations and user requirements makes flexible support for different
95 APIs essential. Of the three current options, the pixel data microservice API has significant
96 advantages from the client's perspective, but not enough OMERO servers have installed the
97 necessary microservice to make it a default. Furthermore, alternatives might become widely
98 adopted in the future, such as the [OMERO Zarr Pixel Buffer](#), and our extension is designed to
99 be able to incorporate these if required.

100 Instructions for the QuPath OMERO extension can be found on a dedicated page of the
101 [QuPath documentation](#). We have also provided javadoc comments for all public fields and
102 methods. The javadocs are installed along with the extension and are available via QuPath's
103 built-in Javadoc viewer to help users write their own scripts.

104 Research impact statement

105 The QuPath OMERO extension was initially released on February 2024. Since then, it has
106 evolved through contributions, bug reporting, and feature requests.

107 The latest extension (v0.2.4) was downloaded more than 4,000 times within the first three
108 months since release. Community enthusiasm for connecting QuPath and OMERO is also
109 reflected in numerous resources and articles that reference the current and previous extensions,
110 while also providing additional documentation ([Biolmaging and Optics Platform, 2026](#); [Burri, 2022](#); [Fan, 2023](#); [Gkotsoulas et al., 2025](#); [Open Microscopy Environment, 2024](#); [Wicks et al., 2023](#)).

113 The usefulness of this work is not limited to users at institutions where images are managed
114 using OMERO. Because OMERO is also used by major public imaging resources, such as the
115 IDR, the extension can help anyone access these resources and explore the data within QuPath.

116 It can further serve as a template for developers to create QuPath extensions that connect to
117 other image management systems.

118 AI usage disclosure

119 No generative AI was used to write the software, documentation or paper.

120 Acknowledgements

121 This work was supported by the Wellcome Trust [223750/Z/21/Z]. This project has been
122 made possible in part by grant number 2021-237595 from the Chan Zuckerberg Initiative DAF,
123 an advised fund of Silicon Valley Community Foundation.

124 We thank Melvin Gelbard and Rémy Dornier for contributions to earlier QuPath OMERO
125 extensions, and the BIOP team (especially Rémy) for extensive testing and feedback on the
126 current work.

127 For the purpose of open access, the author has applied a Creative Commons Attribution (CC
128 BY) licence to any Author Accepted Manuscript version arising from this submission.

129 References

- 130 Allan, C., Burel, J.-M., Moore, J., Blackburn, C., Linkert, M., Loynton, S., MacDonald,
131 D., Moore, W. J., Neves, C., Patterson, A., Porter, M., Tarkowska, A., Loranger, B.,
132 Avondo, J., Lagerstedt, I., Lianas, L., Leo, S., Hands, K., Hay, R. T., ... Swedlow, J.
133 R. (2012). OMERO: Flexible, model-driven data management for experimental biology.
134 *Nature Methods*, 9(3), 245–253. <https://doi.org/10.1038/nmeth.1896>
- 135 Bankhead, P., Loughrey, M. B., Fernández, J. A., Dombrowski, Y., McArt, D. G., Dunne, P.
136 D., McQuaid, S., Gray, R. T., Murray, L. J., Coleman, H. G., James, J. A., Salto-Tellez,
137 M., & Hamilton, P. W. (2017). QuPath: Open source software for digital pathology image
138 analysis. *Scientific Reports*, 7(1), 16878. <https://doi.org/10.1038/s41598-017-17204-5>
- 139 Biolmaging and Optics Platform. (2026). *OMERO - QuPath. BIOP wiki resources*. <https://wiki-biop.epfl.ch/en/data-management/omero/qupath>
- 140
- 141 Burri, O. (2022). QuPath: Taking full control through workflows, scripting, and extensions.
142 *Analytical Science Article DO Series*, 24, 28–31. [https://analyticalscience.wiley.com/](https://analyticalscience.wiley.com/content/article-do/qupath-taking-full-control-through-workflows-scripting-and-extensions)
143 [content/article-do/qupath-taking-full-control-through-workflows-scripting-and-extensions](https://analyticalscience.wiley.com/content/article-do/qupath-taking-full-control-through-workflows-scripting-and-extensions)
- 144 Fan, F. (2023, April 22). *Using OMERO in QuPath to have remote access to your data.*
145 *Tidbits from along the way*. [https://andrewjanowczyk.com/using-omero-in-qupath-to-](https://andrewjanowczyk.com/using-omero-in-qupath-to-have-remote-access-to-your-data/)
146 [have-remote-access-to-your-data/](https://andrewjanowczyk.com/using-omero-in-qupath-to-have-remote-access-to-your-data/)
- 147 Gkotsoulas, D. G., Schönenberger, L., Leupold, J., Callegari, I., Bahn, E., Stadelmann, C.,
148 Elverfeldt, D. von, Kiselev, V. G., Weigel, M., & Granziera, C. (2025). Bridging mesoscopic
149 and microscopic scales in multiple sclerosis: Post mortem brain block multi-contrast 9.4T
150 MRI and histology quantification. *STAR Protocols*, 6(4). [https://doi.org/10.1016/j.xpro.](https://doi.org/10.1016/j.xpro.2025.104203)
151 [2025.104203](https://doi.org/10.1016/j.xpro.2025.104203)
- 152 Open Microscopy Environment. (2024). *QuPath. OMERO guide*. [https://omero-guides.](https://omero-guides.readthedocs.io/en/latest/qupath/docs/index.html)
153 [readthedocs.io/en/latest/qupath/docs/index.html](https://omero-guides.readthedocs.io/en/latest/qupath/docs/index.html)
- 154 Wicks, M. N., Glinka, M., Hill, B., Houghton, D., Sharghi, M., Ferreira, I., Adams, D.,
155 Din, S., Papatheodorou, I., Kirkwood, K., Cheeseman, M., Burger, A., Baldock, R.
156 A., & Arends, M. J. (2023). The comparative pathology workbench: Interactive visual

157 analytics for biomedical data. *Journal of Pathology Informatics*, 14, 100328. <https://doi.org/10.1016/j.jpi.2023.100328>
158
159 Williams, E., Moore, J., Li, S. W., Rustici, G., Tarkowska, A., Chessel, A., Leo, S., Antal,
160 B., Ferguson, R. K., Sarkans, U., Brazma, A., Carazo Salas, R. E., & Swedlow, J. R.
161 (2017). Image data resource: A bioimage data integration and publication platform. *Nature*
162 *Methods*, 14(8), 775–781. <https://doi.org/10.1038/nmeth.4326>

DRAFT